

Technical Specification

| Stapler | Description | Reloads | Color | Open length | Close length | Number of firings | Number of staples | Rows of staples | Staple Line Length |
|---------|------------------------------|------------|-------|-------------|--------------|-------------------|-------------------|-----------------|--------------------|
| MALS30 | Mirus Auto Linear Stapler 30 | MLSR30-3.5 | Blue | 3.5 mm | 1.5 mm | 8 | 11 | 2 | 30 mm |
| | | MLSR30-4.8 | Green | 4.8 mm | 2.0 mm | 8 | 11 | 2 | 30 mm |
| MALS45 | Mirus Auto Linear Stapler 45 | MLSR45-3.5 | Blue | 3.5 mm | 1.5 mm | 8 | 15 | 2 | 45 mm |
| | | MLSR45-4.8 | Green | 4.8 mm | 2.0 mm | 8 | 15 | 2 | 45 mm |
| MALS60 | Mirus Auto Linear Stapler 60 | MLSR60-3.5 | Blue | 3.5 mm | 1.5 mm | 8 | 21 | 2 | 60 mm |
| | | MLSR60-4.8 | Green | 4.8 mm | 2.0 mm | 8 | 21 | 2 | 60 mm |
| MALS90 | Mirus Auto Linear Stapler 90 | MLSR90-3.5 | Blue | 3.5 mm | 1.5 mm | 8 | 33 | 2 | 90 mm |
| | | MLSR90-4.8 | Green | 4.8 mm | 2.0 mm | 8 | 33 | 2 | 90 mm |

Meril's Global Presence

Meril Tibbi Cihazlar

İmalat ve Ticaret A.Ş.
İçerenköy Mah. Çetinkaya Sok.
Prestij Plaza No:28, Kat:4 Ataşehir,
34752 İstanbul / Turkey
T : +90 216 641 44 24
F : +90 216 641 44 25

Meril Medical LLC.

Nauchnyi Proezd 19,
Moscow , Russia – 117 246.
Office - +7 495 772 7643

Meril UK Subsidiary

Crowe UK LLP, St. Brides
House, 10 Salisbury Square,
London EC4Y 8EH, UK

Meril GmbH.

Bornheimer Strasse 135-137,
D-53119 Bonn.
Germany.
T : +49 228 7100 4000
F : +49 228 7100 4001

Meril, inc.

2436 Emrick Boulevard,
Bethlehem, PA - 18020
T : (610)500-2080
F : (610)317-1672

Meril South America

Doc Med LTDA Al. dos Tupiniquins,
1079 - Cep: 04077-003 - Moema.
Sao Paulo. Brazil.
T : +55 11 3624 5935
F : +55 11 3624 5936

EU Representative.

Obelis S.A.Bd, General Wahis 53,
1030, Brussels, Belgium.
T : +32 2 732 5954
F : +32 2 732 6003
E : mail@obelis.net

Meril SA Pty Ltd

102, 104, S101 and S102, Boulevard
West Office Park, 142 Western
Service Road, Erf 813
Woodmead Extension 17 Sandton,
Johannesburg – 2191 South Africa
T : +27 11 465-2049
F : +27 86 471 7941

Meril Endo Surgery Pvt. Ltd

Third Floor E1-E3, Meril Park,
Survey Number 135/2/B & 172/2
Muktanand Marg, Chala,
Vapi 396 191, Gujarat. India.
T : +91 260 240 8000
F : +91 260 305 2125
E : askinfo@merillife.com

Meril China Co. Ltd.

2301b,23f, Lixin Plaza,
no 90, South Hubin Road,
Xiamen, China.
T : 0086-592-5368505
F : 0086-592-5368519

Meril Bangladesh Pvt. Ltd.

Union Heights-01,
55-2, Bir Uttam Qazi
Nuruzzaman Sarak,
West Panthapath,
Dhaka-1205, Bangladesh
T: + 88 02 9122714/24

Meril Australia Pvt. Ltd.

Level 14, 60 Martin Place
Sydney NSW 2000, Australia
T : +61 28248 5800
F : +61 28248 5899

Meril Life Sciences Pvt. Ltd.

301, A - Wing, Business Square,
Chakala, Andheri Kurla Road,
Andheri East, Mumbai 400 093
T : + 91 22 39350700
F : + 91 22 39350777
E : askinfo@merillife.com

**For more information about Endo Surgery products
Please contact your local representative.**

Please see the package inserted for complete device description, product selection information, indications, contraindications, precautions, adverse effects, warning, materials, sterilization and patient guidance associated with the Endo Surgery.

MIRUS™ Linear Stapler



MIRUS™ Linear Stapler




Mirus Linear Stapler (without dial) is engineered to facilitate precise transection & resection in various abdominal & thoracic surgeries.

Mirus Linear Stapler gives two staggered rows of titanium staples in sizes of 30mm, 45mm, 60mm and 90mm; with reload units available for both standard tissue (3.5mm) and thick tissue (4.8mm). The device might be reloaded up to seven times for a maximum of eight firings.

The ergonomic design with single handed closing & firing lever facilitates user comfort and control. Staple formation creates the perfect 'B' shape on each and every firing for improved clinical performance.



Dual Function Handle
For repositioning closure & stapling

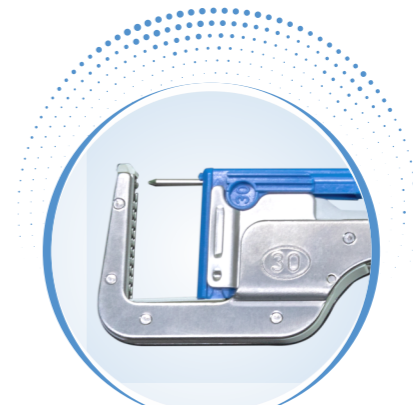

Tissue Retaining Pin
Avoids tissue slip-out



Cartridges
For varied tissue thickness



Cutting-Guide
For safe & precise transection



Quick Release Button
For easy jaw opening



Retaining Pin Knob
Advances tissue retaining pin



Closing & Firing Handle

Handle